


INSTRUCTIONS

Correct Use of the DrySpot biovis® Sample Card




To ensure that your test results are not affected, please observe the following instructions when using the DrySpot sample card:


1. Your healthcare provider should already have attached a barcode to the DrySpot sample card for you. If this is not the case, please label the sample card with your **name and the date of sample collection**.

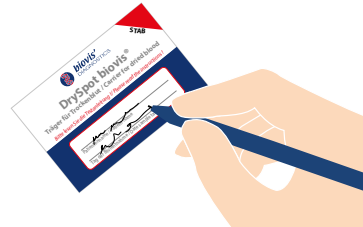


2. Important: Avoid direct sunlight. Place all sampling materials, including the DrySpot sample card, on a clean, dry surface, such as a clean sheet of paper, to prevent contamination that could affect the test results.

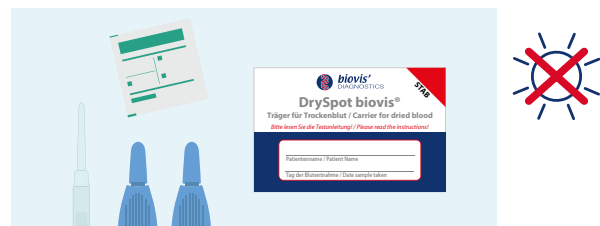


3. The fields for **serotonin and fatty acids** contain different stabilizers and may therefore develop a brownish discoloration – **this is normal**. Always use the appropriately labelled field for the respective test.

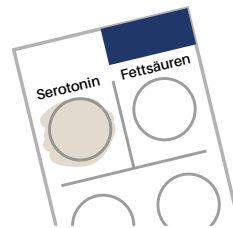


4. Using a Minivette, for example, place the collected blood **in the centre of the pre-punched field**. Please make sure that the pre-punched fields are **completely filled with blood**.

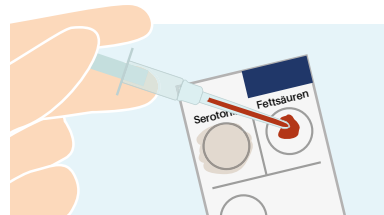

5. Allow the DrySpot sample card to dry for at least **1-2 hours at room temperature in a dark place** (e.g. a darkened room).

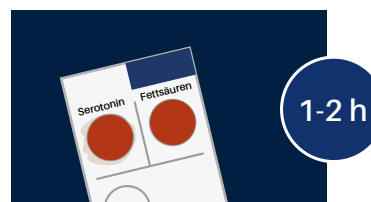

6. Once the samples are completely dry, place the sample card in the **protective sleeve**.

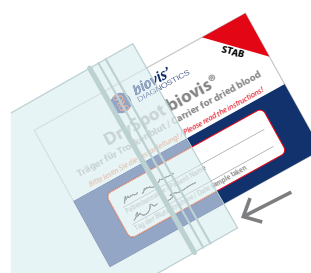

Thank you.